




Bio-Therapeutic

TECHNOLOGY ACCELERATED SKIN CARE™

## bt-micro® quick start guide

1. Read operation manual prior to using the equipment.
2. Do not use the unit if any contraindications are present.

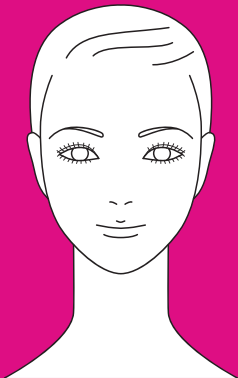


3. Apply appropriate bt-ceuticals® product.

4. Press the **POWER** button for two to three seconds.



5. Press appropriate setting:  
Once for Low, twice for High.  
**EXFOL** to exfoliate  
**SERUM** or **CREAM** to apply serums or creams  
**CLEAR** for clearing.

6. To exfoliate:  
Hold unit in hand, buttons towards palm of hand, applicator end up.  
Glide applicator in a gentle upward motion along the skin, wipe away residue.

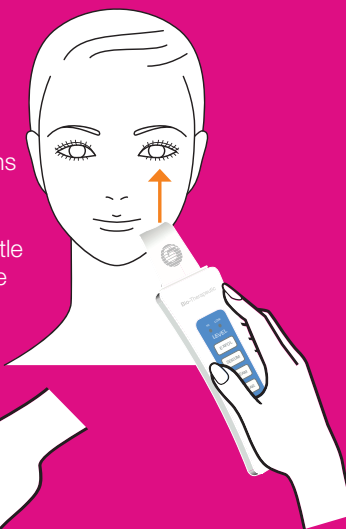

7. To apply DECONGEST, serum, or creams: Hold by the metal sides of unit, buttons away from palm of hand, applicator end down with your free hand 1 - 2 inches away from the applicator to keep skin taut. Glide applicator in an upward, gentle motion.

**NOTE:** Before plugging in your bt-micro mini USB charger, confirm that both the cord and the bt-micro are aligned as shown in the image below with the wider side on the bottom.



mini USB charger



bottom of bt-micro

The bt-micro is a cosmetic exfoliation system. The bt-micro does not penetrate the living layers of the skin, and is not intended in the diagnosis of disease or other conditions, or in the cure, mitigation, treatment, or prevention of disease, and is not intended to affect the structure of the human body. Items may vary without notice and are not available in all countries. Please contact Bio-Therapeutic for updated information.