Osaki

OS-AA15 Neck Massager

IMPORTANT WARNING

- To reduce the risk of electric shock, do not remove the cover. There are no user-serviceable parts inside.
- To reduce the risk of fire or electric shock, do not expose this product to rain or moisture.

The lightning flash with arrow symbol inside an equilateral triangle is intended to alert the user to the presence of uninsulated "hazardous voltage" within the product's enclosure, which may be of sufficient magnitude to induce a risk of electric shock.

The exclamation point inside an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the product.

IMPORTANT: To reduce the risk of electric shock, make sure the plug fits fully into the outlet. If not, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way.

IMPORTANT SAFETY INSTRUCTIONS

When using this massager, basic precautions should always be followed.

A. POTENTIAL DANGERS

To reduce the risk of electric shock:


- Always unplug this massager from the electrical outlet immediately after use and before cleaning.
- Do not reach for a massager that has fallen into water. Unplug it immediately.
- Do not use the massager while bathing or showering.
- Do not place or store the massager where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquids.
- Never use pins or other metallic fasteners with this massager.
- Carefully examine the covering before each use. Dispose of the massager if the covering shows any signs of deterioration and causes blistering and racking on the user's skin.
- Keep the massager dry and do not operate it in a wet or moist environment.

B. WARNINGS

To reduce the risk of burns, fire, electric shocks, or injury to persons:

- The massager should never be left unattended when in use.
- Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shocks, or injury.
- Close supervision is necessary when this massager is used by or near children, dependents, or disabled persons.
- Use this massager only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this massager under the following circumstances: the massager has a damaged cord or plug; the massager has been damaged or encountered a heavy fall; the massager has been dropped into water.
- Do not carry this massager by its power cord or use the cord as a handle.
- Do not crush or fold this massager when storing away.
- Keep the cord away from heated surfaces.
- Never drop or insert any foreign object into any openings.
- Do not use outdoors. This massager is intended for household and indoor use only.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not operate in the presence of explosive or flammable fumes.
- To disconnect, set all controls to the off position and remove the plug from the outlet.
- Do not overload the electrical outlet. Use only the power outlet type as indicated in the "Specifications" section of this manual.
- To avoid the risk of electric shock, do not disassemble or attempt to repair the massager. Incorrect repair can cause risk of electric shock or injuries when the unit is used.
- Never remove the plug from the outlet by pulling the power cord.

Product Description

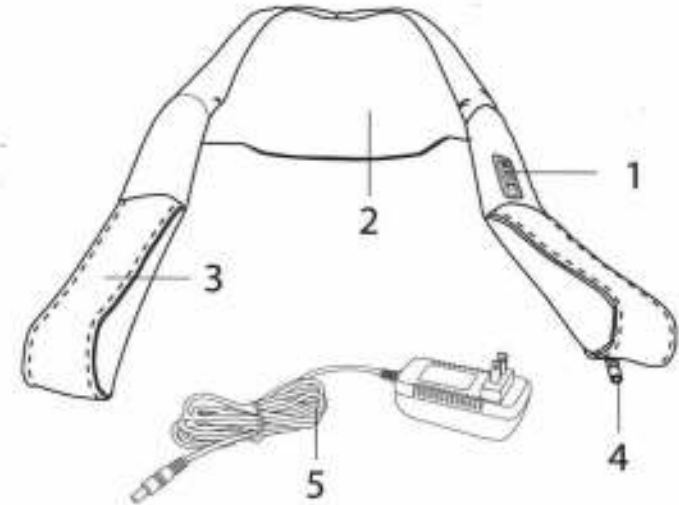- | | |
|------------------------|------------------|
| 1. Control Panel | 4. DC Power Port |
| 2. Shiatsu Nodes | 5. Power Adapter |
| 3. Positioning Handles | |

Operating Instructions

1. Grip the positioning handles and drape the main unit around your neck with the Shiatsu nodes at the center of your back.
2. Press and hold the  ON/OFF Button for 3 seconds and the massager immediately starts in Auto-Reversing Mode. Default Shiatsu intensity will be high.

3. Press the **[P]** Intensity Setting Button to adjust the intensity. The default intensity will be high. Press the Intensity Setting Button to advance the massager to Medium intensity, and press it again for Low intensity.


4. Press the **[C]** Auto Reversing Mode Button to switch the Shiatsu Nodes between Auto Reversing Mode and Single Direction Mode. The default mode is Auto Reversing Mode, which will change mode direction every minute. every one minute.

5. Press the **[W]** Heat Button to activate heat. Press it again to deactivate heat.

6. Use the positioning handles to reposition the Shiatsu nodes on top of the needed treatment area.

7. Press the **[ON/OFF]** button at any time to turn off the massager.

8. The massager will shut off automatically after 15 minutes of continuous use.

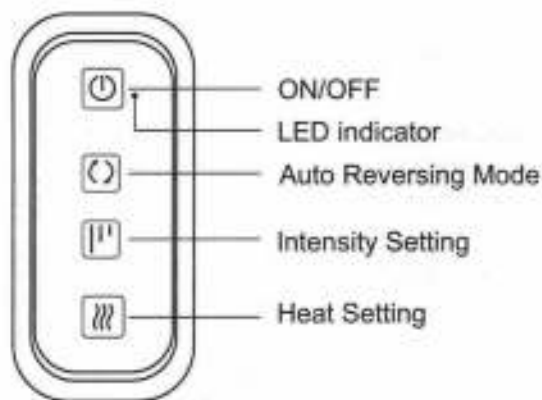Note: It is NOT recommended to use the heat function for more than 15 minutes at the time. Allow the massager to cool down between heat cycles.

Treatment Areas

Maintenance and Care

1. Place the massager in a safe, cool, and dry place when not in use.
2. After unplugging, wipe with a soft, damp cloth using non-abrasive cleaners.
3. Never submerge the massager in any liquids.
4. Keep away from all solvents and harsh detergents.
5. DO NOT attempt to repair this massager. There are no user-serviceable parts.

Specifications

Input: 100V-240V~ 50/60Hz
Rated voltage: DC 12V == 24W

Disposal

Disposal of the Appliance

Under no circumstances should you dispose of the appliance in normal domestic waste.

This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance via an approved disposal company or your municipal waste facility. Please observe currently applicable regulations. Please contact your waste disposal center if you need any further information.

Disposal of Batteries

Used batteries may not be disposed of in household waste.

All consumers are statutorily required to dispose of batteries at a collection point provided by their local municipality or retail store.

The purpose of this obligation is to ensure batteries are disposed of in a non-polluting manner. Only dispose of batteries when they are fully discharged.

The appliance's packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

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IMPORTANT INSTRUCTIONS DE SECURITE

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Safety advice

Read all safety advices and instructions. Non-observance of the safety advices and instructions may cause electric shock, fire and/or bad injuries.

KEEP ALL SAFETY NOTICES AND INSTRUCTIONS FOR FUTURE REFERENCE!

-  **WARNING!**
DANGER TO LIFE AND RISK OF ACCIDENT FOR INFANTS AND CHILDREN! Never leave children unsupervised with the packaging materials. The packaging material poses a suffocation hazard. Children frequently underestimate the dangers. Please always keep the product out of the reach of children.
- ▲ **CAUTION!** This product is not a toy for children! Children are not aware of the dangers associated with handling electrical products.
- This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental

capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved. Children should not clean or perform user maintenance unsupervised.

- ▲ **CAUTION!** The product has a heated surface. Persons insensitive to heat must be careful when using it.
- ▲ **WARNING!** Avoid the hazard of life by electric shock:
 - Do not insert needles or pointed objects into the product.
 - Never immerse the product into water or any liquid.
 - Never touch the product, or the mains power adaptor with wet hands or use it in very humid environments.
 - Only use the provided mains power adaptor as described in Technical Data section.
 - Make sure that the details on the rating label are in conformity with your mains voltage.
 - In the event of malfunctions, remove the mains power adaptor from the mains socket.
 - Do not open the product or attempt to repair it yourself.

Defective products should only be repaired by qualified personnel.

- Always remove the mains power adaptor from the mains socket after each use and before clean the product.
- Keep the power cord away from hot surfaces, sharp edges and mechanical forces.
- Never let the power cord become a trip hazard, that no-one can get tangled in it or tread on it.
- Allow the product to cool down before moving it to a different location.
-  **WARNING! HAZARD OF FIRE!** Do not operate the product close to or below curtains, wall cupboards or other inflammable materials.
- Do not use the product on swollen, burnt, inflamed, diseased or injured regions of the skin and body areas. If in doubt, obtain medical advice before using the product.
- Never use the product if you have injuries or pain in the neck.

● Storage

- Let the product completely cool down before storage.
- Do not place other objects on top of the product.
- If the product is not used for a certain period of time, store the product in a clean and dry location.

● Disposal

The packaging is made entirely of recyclable materials, which you may dispose of at local recycling facilities.

Observe the marking of the packaging materials for waste separation, which are marked with abbreviations (a) and numbers (b) with following meaning: 1-7: plastics / 20-22: paper and fibre-board / 80-98: composite materials.

The product and packaging materials are recyclable, dispose of it separately for better waste treatment. The Triman logo is valid in France only.

Contact your local refuse disposal authority for more details of how to dispose of your worn-out product.

To help protect the environment, please dispose of the product properly when it has reached the end of its useful life and not in the household waste. Information on collection points and their opening hours can be obtained from your local authority.

● Warranty

The product has been manufactured to strict quality guidelines and meticulously examined before delivery.