

Osaki

J1/J5 Jade Mat


Operation Instructions

SAVE THESE INSTRUCTIONS



The Osaki Jade Mat is a healthcare product, which is adopting advanced electronic heating technology, using the healthcare principle of heated jade releasing a variety of elements. This process forms natural magnetic and resonance absorption, the far infrared light improves blood circulation, promotes metabolism, and eliminates waste. This thermal effect, which fundamentally eradicate disease, enhance immunity, realize physiotherapy, and improves healthcare, hairdressing, & fitness.


Safety Warning

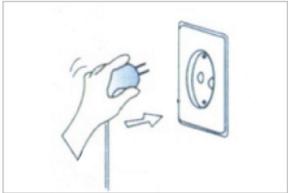

Confirm voltage is 110-120v & 50/60Hz.




Confirm whether power cord is connected to the power socket correctly.

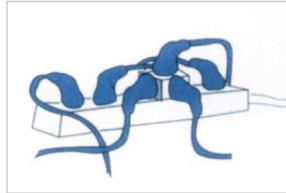

Please don't use it with other electric facilities, this is a hazard.



The perceptually handicapped should only use this product under proper supervision.

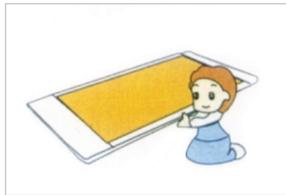


Be careful with the heated surface. Children should be supervised to ensure that they do not play with the appliance.



If there is something wrong with the mat, please stop using immediately and contact dealer.



After using, please turn off the power switch to avoid surges and fire.

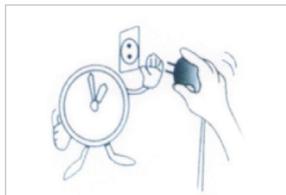


There must be at least a one minute interval when you pull out power cord and connect it again.



When you clean the mattress, please disconnect from power first, as to avoid getting an electric shock.

Name and function components

- 1st Layer:** High-grade leather
- 2nd Layer:** Dense mesh cloth
- 3rd Layer:** Jade, germanite or tourmaline stone
- 4th Layer:** Bamboo carbon fiber cotton
- 5th Layer:** Glass fiber
- 6th Layer:** Heating Layer
- 7th Layer:** Heat-sensitive device layer
- 8th Layer:** Ecological nonwoven cloth layer
- 9th Layer:** Heat gathering and releasing layer
- 10th Layer:** Ecological compression cotton
- 11th Layer:** Thermal insulation
- 12th Layer:** The waterproof eleusine cloth

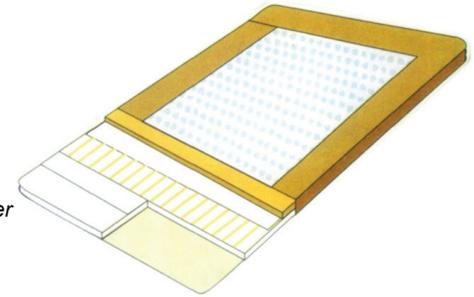

1st Layer

This leather layer is elegant, durable and easy to clean.

2nd Layer

This layer uses high density mesh cloth with elegant embroidery to hold the jade, germanite and tourmaline stones.

3rd Layer

Jade, germanite and tourmaline stone can release infrared light and negative ions, which can promote blood circulation and improve metabolism.

4th Layer

1. The absorption capacity of bamboo carbon is more than charcoal, it can absorb dust, decompose smell and deodorize
2. Releasing infrared light, and storing the heat energy
3. The bamboo carton can adjust humidity
4. The emission concentration of anion is very high, which is beneficial to health
5. High mineral content

5th Layer

This material can cut off electromagnetic wave.

6th Layer

This heat layer uses standard heating material, after heating, the temperature can rapidly rise. This heating material does not easily age, and is resistant to high temperatures and folding.

7th Layer

This heat-sensitive device can prevent inside from overheating, once the temperature is too high, this device will automatically power off.

8th Layer

This layer fixes the heating lines.

9th Layer

This layer stores the heating energy of far infrared light, so that can balance heating area and continually release heat.

10th Layer

Ecological compression cotton layer.

11th Layer


Stops conducting heat to bottom, and prevents overheating.

12th Layer


Durable, prevents wear & tear, resists tension and protects the mattress.

Name and function components


J5 Mat Remote Control

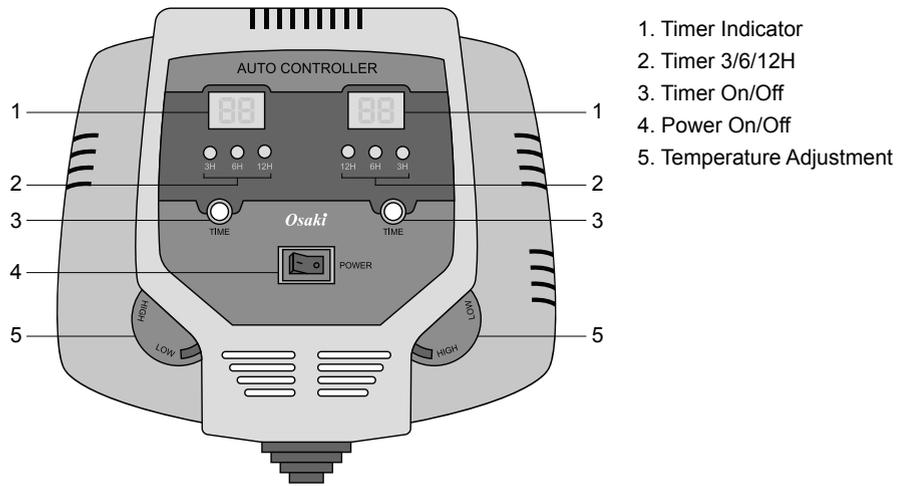

Remote Control Operation Guide

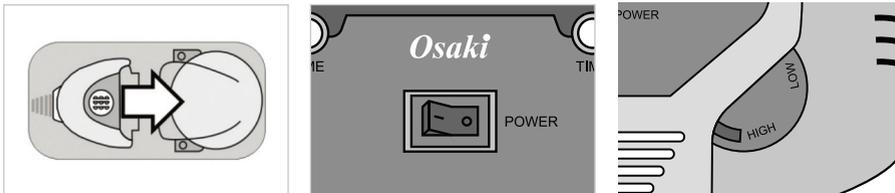


1. Connect the remote controller to the mattress.
2. Turn on the power.
3. Adjust temperature of the mattress

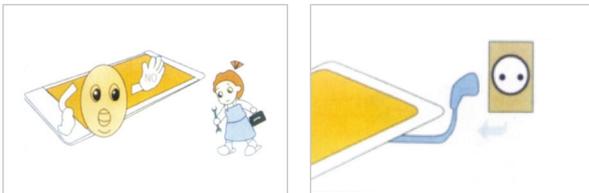


4. When you use this product, the remote controller should always be put outside the mat.
5. After using, turn off the power switch and unplug the power cord.

Name and function components


J1 Mat Remote Control

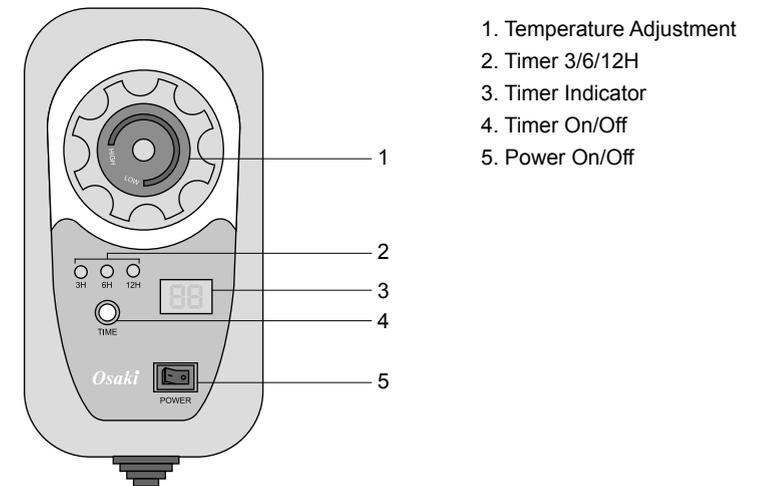

Remote Control Operation Guide

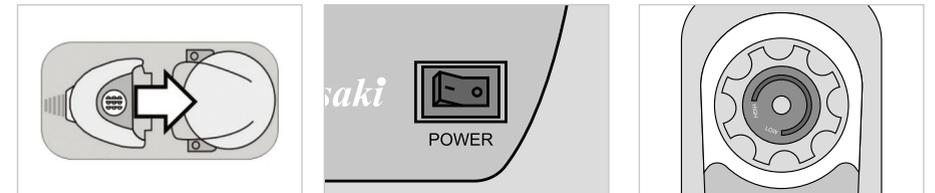

1. Connect the remote controller to the mattress.
2. Turn on the power.
3. Adjust temperature of the mattress

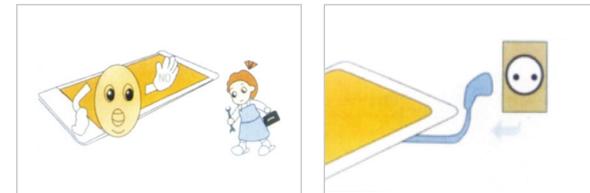

4. When you use this product, the remote controller should always be put outside the mat.
5. After using, turn off the power switch and unplug the power cord.

Method of Usage

During Use

1. After using, turn off the power switch and pull out power cord.
2. Avoid contact with the substances which can be easily set on fire.
3. Avoid to putting product in direct sunlight, liquids and impurities should never enter the interior of product.
4. Product should be stored in dry and well-ventilated place (10°C-25°C).
5. Don't force power cord into socket, this can lead to a total malfunction.
6. Please don't overload the product.

After Use

1. You should clean the mattress thoroughly after use.
2. Mattress should be heated in high temperature for a little while, to evaporate the humidity of the mattress.
3. Folding mattress: the surface of mattress should be inward-folded three times (the part of power socket connection should be folded inward first) for storage.
4. The mattress should be stored in the dustproof bag.
5. The mattress should be stored in dry and well-ventilated place.
6. After storing for a long period of time, you should take it out and check the product's condition.
7. Don't stack heavy things on the mattress when you store it.

Cleaning Method

1. When the mattress needs cleaning and maintenance, please turn off power first.
 2. The surface of the product should be cleaned by dry cloth.
 3. Use a cleaner to quickly absorb the dust between the stones.
 4. Do not use a volatile solvent to clean mattress.
 5. Keep product dry, avoid to contact with water or other liquids.
 6. When you finish cleaning, let mattress air dry.
 7. Please use dry washcloth to clean the surface of controller.
- *To guarantee the life of the product, you should pay attention close attention to these details.

Method of Usage

About Jade Stone

Jade is rare gem. It includes more than 30 kinds of micro elements, just like selenium, Zinc, Nickel, Copper, Manganese, Cobalt, Magnesium, Calcium and so on. These micro elements are beneficial for human health.

The waves of these elements are similar to the human cells. With the jade produces the wave, and the human cells produce a resonance that makes the cell tissue healthier, promotes blood circulation, promptly eliminate body waste. The magnesium is one of the major components of Chlorophyll, magnesium also exists in human cells. Jade can make drinking water and natural water into active water, which can help the body to enhance immunity.

Jade can produce special photoelectric effect, focus on energy, and forms a "field", which produce resonance with the human body, thus making the body function more coordinated.

About Germanite Stone

Germanite stones include large amounts of organic germanium, precious jade powder, pearl powder, Websterite and medicine materials. These minerals are good for human health. Germanite are crystallized under 1350°C. Through heated carbon fiber and releasing a large amount of anion, the stones achieve significant effects of traditional Chinese medicine magnets.

About Tourmaline Stone

Tourmaline has more than 10 kinds of microelements, just like magnesium, iron, aluminum, boron and so on. These micro elements are good for human health. It is a kind of special structure of acute crystals, which can produce electrical ion and permanently release negative ions and infrared light. The biggest characteristic of Tourmaline stone is to produce a permanent 0.06A electric current, it can improve blood circulation.

Product Specification

Product Name: J1/J5 Jade Mat

Rated Voltage: AC 110 - 120V

Rate Power: 440W

Rated Frequency: 50/60Hz

G.W / N.W: 69/7lbs / 64.2 lbs