



# WATER FLOSSER ADVANCED

## QUICK START GUIDE



### ON-DEMAND BUTTON

Activates water  
stream only when  
button is pressed  
down for increased  
control and  
less mess

### ON/OFF BUTTON FOR CONTINUOUS STREAM

Activates water  
stream until button is  
pressed again to stop  
stream


intense

medium

sensitive

MODE  
SELECTION  
Selects preferred  
stream intensity

# GET STARTED

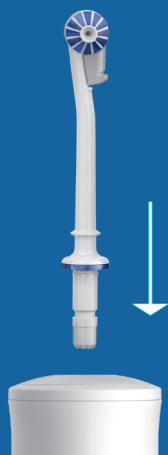

1

PUSH UNTIL BASE OF NOZZLE TOUCHES HANDLE.

2

CHARGE UNTIL FLASHING LIGHT TURNS OFF.



3

FILL WATER RESERVOIR WITH LUKEWARM WATER.

4

AFTER USE: OPEN TANK & DISCARD RESIDUAL WATER. LEAVE TANK DOOR OPEN TO DRY. RESERVOIR CAN BE REMOVED AND IS DISHWASHER SAFE.\*



\*ONLY RESERVOIR IS DISHWASHER SAFE, NOT THE HANDLE ITSELF.

## AQUAFLOSS NOZZLE




PUSH SWITCH UP  
FOR TARGETED  
STREAM.




PULL SWITCH  
DOWN FOR  
ROTATIONAL  
STREAM.

## PRECISION JET NOZZLE

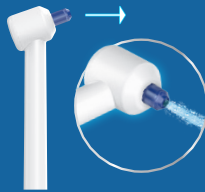

CLICK OUT  
FOR SINGLE  
HIGH-PRESSURE  
STREAM.

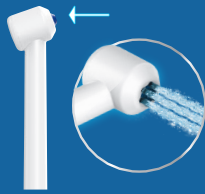

CLICK IN FOR  
MULTI-JET  
STREAM.

## FEATURES



### TWO WAYS TO START



PRESS TO START A CONTINUOUS STREAM



OR



PRESS TO START STREAM, RELEASE TO STOP STREAM.



INTENSITY SELECTION

- INTENSE
- MEDIUM
- SENSITIVE

# HOW TO USE

1. SELECT NOZZLE AND CHOOSE STREAM SETTING.
2. PRESS MODE SELECTION BUTTON TO CHOOSE DESIRED INTENSITY.
3. PUT NOZZLE IN MOUTH. LEAN OVER SINK, CLOSE LIPS TO PREVENT SPLASHING.
4. PRESS THE POWER BUTTON OR THE ON-DEMAND BUTTON.
5. KEEP UNIT UPRIGHT DURING USE FOR CONTINUOUS WATER FLOW.
6. AIM STREAM AND FOLLOW GUMLINE AT 90°. PAUSE BRIEFLY BETWEEN TEETH.



AQUAFLOSS  
NOZZLE



PRECISION JET  
NOZZLE



ORDER MORE NOZZLES ONLINE

LEARN MORE AT [ORALB.COM](https://www.oralb.com)

USE BEFORE BRUSHING FOR BEST RESULTS